

File 351:Derwent WPI 1963-2004/UD,UM &UP=200474

(c) 2004 Thomson Derwent

\*File 351: For more current information, include File 331 in your search.

Enter HELP NEWS 331 for details.

Set	Items	Description
---	-----	-----
? s	pn=ep	591985
	S1	1 PN=EP 591985
? t1/7		

1/7/1

DIALOG(R)File 351:Derwent WPI

(c) 2004 Thomson Derwent. All rts. reserv.

009840283 \*\*Image available\*\*

WPI Acc No: 1994-120139/ 199415

**Glenoid alignment guide for aligning drill bit - has retractor plate to displace muscles, drill guide element which moves along plate, and curved finger to engage glenoid**

Patent Assignee: MINNESOTA MINING & MFG CO (MINN )

Inventor: SHEARER J R; SHELLEY P

Number of Countries: 009 Number of Patents: 003

Patent Family:

Patent No	Kind	Date	Applicat No	Kind	Date	Week
EP 591985	A1	19940413	EP 93116322	A	19931008	199415 B
CA 2107726	A	19940410	CA 2107726	A	19931005	199426
US 5437677	A	19950801	US 93128844	A	19930929	199536

Priority Applications (No Type Date): GB 9221257 A 19921009

Cited Patents: 1.Jnl.Ref; US 5030219

Patent Details:

Patent No	Kind	Lan	Pg	Main IPC	Filing Notes
-----------	------	-----	----	----------	--------------

EP 591985	A1	E	10	A61B-017/16	
-----------	----	---	----	-------------	--

Designated States (Regional): CH DE FR GB IT LI SE

US 5437677	A	8	A61B-017/56
------------	---	---	-------------

CA 2107726	A		A61B-017/56
------------	---	--	-------------

Abstract (Basic): EP 591985 A

The guide has a retractor plate (32) with a curved finger (34), narrower than the plate, at its forward end and this finger is for insertion behind the glenoid (12). The tip (10) of the finger engages the bone of the glenoid neck (10) to form a fulcrum. Slidably carried on the retractor plate is a drill guide (40).

A guide device (60) guides a drill bit (62) along the path to make a hole (63) for the glenoid face (26) into the glenoid neck (10). When the tip of the plate engages the bone, the axial path of the drill guide is aligned relative to the neck.

ADVANTAGE - The arrangement ensures accurate placement of implant.

Dwg.3/7

Abstract (Equivalent): US 5437677 A

A glenoid alignment guide for aligning a drill bit relative to the glenoid and the glenoid neck of a patient, comprises a retractor plate for displacing the posterior and superior aspect of the deltoid muscle, and a drill guide slidably movable along the retractor plate. the tip of a narrow, curved finger on the retractor plate is adapted to engage the bone of the glenoid neck to form a fulcrum for levering the muscle and other tissue clear of the glenoid.

When the tip of the retractor plate engages the glenoid neck and the drill guide is advanced along the retractor plate against the

glenoid face, the drill guide is adapted to guide a drill bit to make a hole from the glenoid face into the glenoid neck along the axis of the glenoid neck. This hole is then used to receive a fastener or locating peg of a glenoid component of a shoulder implant.

ADVANTAGE - The glenoid alignment guide is adapted among the other things for improved viewing of the glenoid face. For example, the drill guide may have a periphery and a number of cut away portions around the periphery to improve the view of the glenoid face as the drill guide is advanced. This allows the glenoid face to be viewed through the cut away portion of the drill guide to facilitate centering the drill guide relative to the glenoid face.

Dwg.3/6

Derwent Class: P31

International Patent Class (Main): A61B-017/16; A61B-017/56

?

(19)



Europäisches Patentamt  
European Patent Office  
Office européen des brevets



(11) Publication number:

**0 591 985 A1**

(12)

**EUROPEAN PATENT APPLICATION**

(21) Application number: 93116322.4

(51) Int. Cl.<sup>5</sup>: **A61B 17/16**

(22) Date of filing: 08.10.93

(30) Priority: 09.10.92 GB 9221257


(43) Date of publication of application:  
13.04.94 Bulletin 94/15(84) Designated Contracting States:  
CH DE FR GB IT LI SE(71) Applicant: **MINNESOTA MINING AND  
MANUFACTURING COMPANY**  
3M Center,  
P.O. Box 33427  
St. Paul, Minnesota 55133-3427(US)**Road****P.O. Box 33427****St. Paul****Minnesota 55133-3427(US)****Inventor: Shelley, Philip c/o 2501 Hudson  
Road****P.O. Box 33427****St. Paul****Minnesota 55133-3427(US)**

(72) Inventor: Shearer, John R. c/o 2501 Hudson

(74) Representative: **VOSSIUS & PARTNER**  
**Postfach 86 07 67**  
**D-81634 München (DE)**(54) **Glenoid alignment guide and method of use.**

(57) A glenoid alignment guide (30) for aligning a drill bit relative to the glenoid and the glenoid neck of a patient. The alignment guide comprises a retractor plate (32) for displacing the posterior and superior aspect of the deltoid muscle, and a drill guide (40) slidably movable along the retractor plate. The tip (36) of a narrow, curved finger (34) on the retractor plate is adapted to engage the bone of the glenoid neck to form a fulcrum for levering the muscle and other tissue clear of the glenoid. When the tip of the

retractor plate engages the glenoid neck and the drill guide is advanced along the retractor plate against the glenoid face, the drill guide is adapted to guide a drill bit to make a hole from the glenoid face into the glenoid neck along the axis of the glenoid neck. This hole is then used to receive a fastener or locating peg of a glenoid component of a shoulder implant. The method of using the alignment guide is also disclosed.

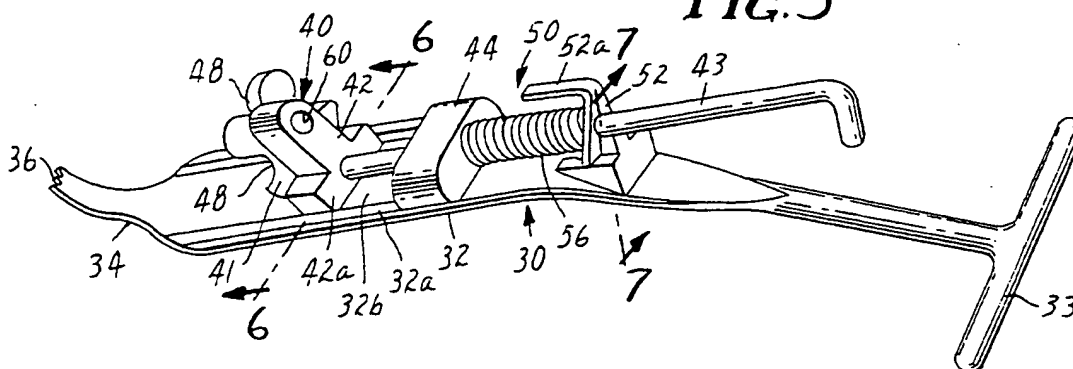**FIG. 3****EP 0 591 985 A1**

This invention relates to a glenoid alignment guide to assist in the preparation of the glenoid during surgery to replace a shoulder joint with an implant.

The preparation of the glenoid during shoulder surgery to receive a plastic or metal implant relies upon the surgeon being able to identify the axis of the bone and locate the implant within the maximum bone stock. Access to the glenoid is limited by the depth of the structure in the wound and the surrounding tissue as well as the surgical technique used.

The anatomy of the glenoid and scapula is such that there is very little bone into which an implant can be placed. The articulating surface of the joint consists of a shallow dishd cartilaginous area bonded by soft tissue of the rotator cuff which stabilizes the humeral head against the glenoid.

The glenoid is roughly oval to pear-shaped and most replacement implants follow this contouring for optimal fit. Whilst the lateral aspect of the glenoid presents a broad surface for articulation, the margins taper rapidly medially to form a narrow neck of about 10 to 15mm thickness from which emerge the coracoid (superiorly) and acromial (posteriorly) processes. The thinning of the glenoid progresses into the scapula bone where the thickness may be as little as 2mm to 3mm, although there is a broader spine running along the inferior margin of the scapula.

The axis of the visible face of the glenoid is easily identified, but does not correspond to the axis of the glenoid neck which is oriented posteriorly behind the glenoid face. The amount of bone in the neck in line with the axis of the visible glenoid is therefore greatly reduced and so it is along the actual axis of the glenoid neck that the anchoring means of the implant needs to be positioned. In addition, disease or degeneration can severely restrict the size of the glenoid, making implant placement more problematic. Further the scapula is a floating bone which is highly mobile, presenting little resistance to pressure, particularly under anaesthesia.

There are several problems for the accurate and satisfactory placement of an implant including:

1. access to the glenoid;
2. identification of the axis of the glenoid;
3. the alignment of the implant against the glenoid face,
4. the alignment of the anchoring means for the implant in the glenoid bone, and
5. the risk of glenoid neck fracture or exposure and penetration of the anchoring means through the wall of the bone.

To date the placement of an implant has depended upon the surgeon visually assessing the axis of the glenoid. He can gain access to and see

the glenoid. However, the deltoid muscle, tendons and other tissue obscure the glenoid neck and it is undesirable to disturb this region unnecessarily. Therefore, the surgeon has had little or no actual knowledge of the anatomy behind the face of the glenoid, particularly the thickness and orientation of the neck which as noted above does not correspond with the axis of the glenoid. There is therefore always a risk that the anchoring device for the implant will either impinge on the sloping back wall of the neck which will lead to unsatisfactory placement of the implant against the glenoid face, or will penetrate the cortical bone which may result in bone fracture or interference with or damage to soft tissue.

The invention therefore has been made with these factors in mind. The invention provides a glenoid alignment guide adapted to retract soft tissue and muscle away for the glenoid, and guide a drill bit into the glenoid in proper alignment with the glenoid neck and glenoid face of a patient.

Generally, the glenoid alignment guide comprises a retractor plate for displacing the posterior and superior aspect of the deltoid muscle, and a drill guide slidably carried on the retractor plate. The retractor plate has a curved finger narrower than the plate defining the forward end of the retractor plate. The finger is adapted for insertion behind the glenoid. The tip of the finger is arranged to engage the bone of the glenoid neck to form a fulcrum for levering the muscle and other tissue clear of the glenoid. The drill guide defines an axial path along which a drill bit is guided. The drill guide is movable along the retractor plate to a forward position for engagement with the glenoid face. Guide means (e.g., a guide hole) is provided in the drill guide for guiding a drill bit along the axial path to make a hole from the glenoid face into the glenoid neck along the axis of the glenoid neck to receive fastening or locating means of a glenoid component of a shoulder implant. The tip of the retractor plate when engaged against the bone of the glenoid neck aligns the axial path defined by the drill guide relative to the glenoid neck.

Preferably, the tip of the finger has serrations adapted for non-slip engagement with the neck of the glenoid, and the tip of the finger is spaced radially from the center of the axial path defined by the guide means by 1.5mm to 2.5mm plus the radius of the drill bit to be guided by the guide means. Most preferably, the tip of the finger is spaced radially from the center of the axial path of the guide means by 2.0mm to 2.5mm plus the radius of the drill bit to be guided by the guide means.

Also, preferably, the glenoid alignment guide is adapted among other things for improved viewing of the glenoid face. For example, the drill guide

may have a periphery and a plurality of cut away portions around the periphery to improve the view of the glenoid face as the drill guide is advanced. This allows the glenoid face to be viewed through the cut away portions of the drill guide to facilitate centering the drill guide relative to the glenoid face.

Most preferably, the drill guide has a generally convex-curved frontal surface adapted to correspond to the generally concave-curved surface of the glenoid face of the glenoid.

In one preferred aspect, locking means is provided for releasably locking the drill guide against the glenoid face of the glenoid. For example, the locking means may comprise a rod projecting rearwardly from the drill guide, a pivotable locking lever having a hole with an elongate cross section through which the rod passes, and resilient urging means for urging the locking lever such that the rod is gripped by the edges of the hole. The arrangement is such that retraction of the drill guide is resisted by the gripping of the rod by the locking lever but is possible if the locking lever is manually pivoted against the urging force of the resilient urging means to release the rod.

Most preferably, the retractor plate has a track mounting the drill guide on the retractor plate for sliding movement of the drill guide along the track in the axial direction. The track constrains the drill guide from rotation relative to the retractor plate, with the drill guide being mounted on the rod to move axially together with the rod to allow advancing the drill guide along the track of the retractor plate by manually moving the rod. For example, the retractor plate has a main body having opposite side edges, and ribs extending from the main body along the opposite side edges to define the track.

In yet another aspect, the drill guide has a rear surface, and a collar is fitted around the drill bit for engaging the rear surface of the drill guide to limit the penetration of the drill bit into the glenoid bone.

Such a guide can assist the surgeon in both obtaining good exposure of the glenoid and location for drilling a hole for fixing an implant. Thus the finger can be relatively narrow and so damage to nerves and blood vessels behind the glenoid can be kept to a minimum as it is inserted and its tip, when engaged with the bone, will identify the thinnest part of the glenoid neck so aligning the drill guide for optimum drilling. Thus the neck of the glenoid is of reasonably consistent width in all patients and once the tip of the finger engages it, one can be sure that the drill guide will be well positioned and ensure that, when the surgeon drills a hole, it will pass along the bone in the neck and not penetrate through the neck.

In a further aspect of the invention, referring to the method of using the glenoid alignment guide, the retractor plate is inserted between the patient's

soft tissue and muscle and the glenoid, and the tip of the finger is placed against the glenoid neck. With the tip of the finger against the glenoid neck, the tissue and muscle is levered aside to expose the glenoid face, by means of the retractor plate. Then, with the tip of the finger remaining against the glenoid neck, the drill guide is advanced to a position where the drill guide abuts the glenoid face. With the tip of the finger remaining against the glenoid neck and the drill guide abutting the glenoid face, the drill guide is used to align a drill bit for drilling a hole along the axial center of the glenoid neck to receive fixing means of the implant.

The invention will now be described, by way of example, with reference to the accompanying drawings, in which corresponding parts are indicated by corresponding reference characters, and in which:

Figures 1 and 2 are views showing a human scapula;

Figure 3 is a perspective view of a glenoid alignment guide according to the invention;

Figure 4 is a side view of the glenoid alignment guide shown in position for aligning the drilling of a hole to receive an implant in the glenoid;

Figure 5 is an enlarged detail of the glenoid alignment guide and glenoid shown in Figure 4 and additionally showing the drilling of the hole;

Figure 6 is an enlarged cross-sectional detail taken along the line 6-6 of Figure 3; and

Figure 7 is an enlarged cross-sectional detail taken along the line 7-7 of Figure 3.

As can be seen from Figures 1 and 2, the neck 10 of the glenoid 12 of a human scapula 14 is very narrow, particularly as viewed in the superior orientation shown in Figure 2. It is therefore important when making a drilled hole 63 into the glenoid 12 to receive the fixing screw of an implant, to ensure that the axis of the drilled hole 63 passes along the actual axis 22 of the glenoid 12, and not along the apparent axis 24 which is the one visible to a surgeon who can only see the articulating face 26 of the glenoid, the rest of the glenoid being covered by the deltoid muscle, tendons and other tissue, etc.

In accordance with the invention a glenoid alignment guide 30 as shown in Figures 3 to 5 is used by a surgeon to assist in identifying this axis 22. This alignment guide 30 includes a retractor plate 32 which is used to displace the posterior and superior aspects of the deltoid muscle to assist in placing the glenoid alignment guide against the glenoid articulating face 26.

At its rear, the retractor plate 32 has a handle 33 to assist in manoeuvring the glenoid alignment guide 30, whilst at its front it has an integral curved finger 34 with a tip 36. This curved finger 34 acts as a fulcrum for the retractor plate 32 when the

surgeon inserts the finger 34 behind the glenoid 12 and can then use the retractor plate 32 to lever muscle and the tissue clear of the glenoid face 26.

The finger 34 is narrower than the retractor plate 34 so as to reduce the damage to nerves and blood vessels which run behind the glenoid face 26. The tip 36 is not, however, pointed, or so narrow that it will penetrate into the bone as the powerful deltoid muscle is levered out of the way by the alignment guide 30. The tip 36 preferably has serrations as shown, or forked (not shown), to provide a secure non-slip engagement with the glenoid neck 10 so that the tip 36 does not slip as the muscle, etc. is being levered clear. The plate 32, therefore, assists the surgeons by ensuring that the glenoid 12 and particularly its face 26 are fully exposed.

The tip 36 of the finger 34 also serves as a guide. The surgeon locates the narrowest part of the glenoid neck with this finger 34. By ensuring that the tip 36 is kept in touch with this part of the bone, he can then ensure that the glenoid alignment guide 30 is located so that an eventual hole 63 drilled using the alignment guide 30, as will be described, will be along the axis 22 (Figure 2), and so make best use of the available bone and will not break through the bone.

A drill guide 40 is slidably mounted on the retractor plate 32. The drill guide 40 includes a front portion 41 and an integral rear portion 42. The retractor plate 32 has upstanding ribs or edges 32a which, as best shown in Figure 6, are slightly above the central surface 32b of the surface of the retractor plate 32 to define a track for slidably mounting the drill guide 40. These edges 32a engage side wing portions 42a extending from the rear portion 42 of the drill guide 40. Attached to and extending rearwardly from the rear portion 42 of the drill guide 40 is a rod or shaft 43 which is slidably mounted in a block 44 fixed to the retractor plate 32. The mounting of the shaft 43 and the engagement of the sides of the drill guide 40 with the edges 32a, ensures that the drill guide 40 is constrained to move axially along the retractor plate 32.

The front face 46 of the drill guide 40 is convex-curved and designed to engage and approximately match the generally concave-curved articulating face 26 of the glenoid. In addition three cut-outs 48 are provided to improve the surgeon's view of the face 26 when the drill guide 40 is brought to abut the face 26 as will be described. Thus, the cut outs 48 allow the surgeon to see the glenoid face and assess the correct size of implant to choose to be fitted. Further, the cut outs are positioned so as to coincide with the positioning of the small stabilizing screws which are used to prevent rotation of the fitted implant, and so facilitate the surgeon's

assessment of whether one or more such screws can be satisfactorily positioned.

In order to lock the drill guide 40 in a selected position, a locking mechanism 50 is provided. This comprises a locking lever or clamp 52 pivotally mounted on the retractor plate 32, having a forwardly-bent portion 52a, and having a hole 54 having an elongate cross section (Figure 7) through which the shaft 43 extends, and a coil spring 56 positioned around the shaft 43 between the block 44 and the clamp 52. The arrangement is such that the coil spring 56 normally pushes the locking clamp 52 in a rearward direction, the edges of the hole 54 clamping the shaft 43 and preventing its rearward movement. The drill guide 40 and shaft 43 can however be pushed forwardly since such movement will release the clamping action of the clamp 52 and temporarily overcome the spring force. In the event that the drill guide 40 is to be retracted, then one presses down on the forwardly-bent portion 52a which moves the locking clamp 52 forwardly against the action of the spring to release the shaft 43, so that the shaft and the drill guide 40 can be withdrawn.

Once the surgeon has positioned the tip behind the glenoid, then using the shaft 43, he pushes the drill guide 40 forwardly until its front face 46 contacts and presses firmly against the articulating face 26 of the glenoid 12, the cut outs 48 allowing the surgeon to see parts of the glenoid surface 26 to help ensure that tissue is not trapped between the face 46 and surface 26. This position is shown in Figures 4 and 5. The locking mechanism 52 then holds the drill guide 40 in that position.

As best shown in Figures 3 and 5, the drill guide 40 has a central through hole 60, and this is used to align a drill bit 62. The through hole 60 has a cross section sized to closely receive the drill bit 62 to provide adequate guidance of the drill bit 62. Therefore the surgeon can now drill a drilled hole 63 in the glenoid as shown in Figure 5.

As can be seen from that figure, the axial path or axis 60a of the through hole 60 along which the drill bit 62 will be advanced is spaced approximately 2mm to 2.5mm plus the radius of the drill bit 62 from that tip 36 of the finger 34. In other words, the spacing 64 (Figure 5) between path of the drill bit 62 and the tip 36 is preferably 2mm to 2.5mm. This ensures that hole 63 is made approximately in the center, in the sense viewed in Figure 5, of the glenoid neck 10.

Since the through hole 60 is sized to closely receive the drill bit 62, one alternative way of expressing this is that the tip 36 of the finger 34 is spaced from the extended central axis 60a of the through hole 60 by a distance approximately equal to the radius of the through hole 60 plus 2mm to 2.5mm.

More specifically, the glenoid neck 10 has a dimension, in the direction of the arrow 65 (Figure 5), which is typically 10mm and almost never more than 12mm. Assuming therefore that the drill bit 62 has a diameter of 5mm, this leaves a thickness of 5mm, or exceptionally, 7mm of bone, and so 2.5mm, i.e., half the typical 5mm of available bone, is the preferred maximum spacing of the tip 36 of the finger 34 from what will be the hole 63 which is drilled. In case the dimension of the glenoid neck 10 may be slightly less than 10mm, however, this spacing 64 should be less than that maximum, e.g., 2mm, and desirably no less than 1.5mm. By choosing the spacing 64, to be 1.5mm to 2.5mm and most preferably approximately 2mm, one can therefore ensure that the drilled hole 63 is approximately centered in the glenoid neck 10 in the sense of the arrow 65.

The alignment of the drilled hole 63 in the sense transverse to the direction of the arrow 64, is not as critical and so small variations of the alignment of the drilled hole 63 in the sense transverse to the plane of the Figure 5, is not usually too critical.

A cylindrical sheath or collar 66 is fitted over the drill bit 12 before it is inserted through the through hole 60 of the drill guide 40. This collar 66 prevents excessive drilling of the glenoid 12 since it will engage the rear face of the drill guide 40 when the depth of the drilled hole 63 is sufficient, normally 20 to 25mm. Thus the drill bit 12 should not be allowed to penetrate through the back of the scapular.

Therefore the through hole 60 of the drill guide 40 guides and aligns the axis of the drilled hole 63 which is drilled and the sheath 66 abuts the rear face of the rear portion of the drill guide 40 when the hole 63 has been drilled to sufficient depth, so limiting the depth of the drilled hole 63. Both arrangements therefore prevent the drill bit 12 and drilled hole 63 from breaking through the neck 10.

The drill guide 40 can thereafter be retracted and the alignment guide 30 removed to allow subsequent steps in the fitting of the implant which may include preparation of the glenoid articulating surface, tapping of the hole 63, and securing of the implant by screwing into the tapped hole 63. These steps form no part of the present invention and so need no further explanation. They are, however, steps which are known and understood by the surgeon.

Co-assigned British Patent Application No. 2 251 795 discloses a preferred shoulder implant, and is hereby incorporated by reference. See, also, the U.S. counterpart to that application, Serial No. 08/035,906, filed March 23, 1993, by John R. Shearer and Philip Shelley, on "Orthopedic Implant", which is a continuation of Serial No.

07/822,619, filed January 17, 1992, which is also hereby incorporated by reference.

As various changes could be made in the above constructions and methods without departing from the scope of the invention, it is intended that all matter contained in the above description or shown in the accompanying drawings be interpreted as illustrative and not in a limiting sense.


## Claims

1. A glenoid alignment guide (30) for aligning a drill bit (62) relative to the glenoid (12) and the glenoid neck (10) of a patient, the glenoid alignment guide (30) comprising:
  - a retractor plate (32) for displacing the posterior and superior aspect of the deltoid muscle, the retractor plate (32) having a curved finger (34) narrower than the plate (32) defining the forward end of the retractor plate (32), the finger (34) being adapted for insertion behind the glenoid (12), the tip (36) of the finger (34) being arranged to engage the bone of the glenoid neck (10) to form a fulcrum for levering the muscle and other tissue clear of the glenoid (10); and
  - a drill guide (40) slidably carried on the retractor plate (32) and defining an axial path along which a drill bit (62) is guided, the drill guide (40) being movable along the retractor plate (32) to a forward position for engagement with the glenoid face (26), and having guide means (60) for guiding a drill bit (62) along the axial path to make a hole (63) from the glenoid face (26) into the glenoid neck (10) along the axis of the glenoid neck (10) to receive fastening or locating means of a glenoid component of a shoulder implant, the tip (36) of the retractor plate (32) when engaged against the bone of the glenoid neck (10) aligning the axial path defined by the drill guide (40) relative to the glenoid neck (10).
2. A glenoid alignment guide (30) according to claim 1 in which the tip (36) of the finger (34) has serrations adapted for non-slip engagement with the neck (10) of the glenoid (12).
3. A glenoid alignment guide (30) according to claims 1 or 2 in which the tip (36) of the finger (34) is spaced radially from the center (60a) of the axial path defined by the guide means (60) by 1.5mm to 2.5mm plus the radius of the drill bit (62) to be guided by the guide means (60).
4. A glenoid alignment guide (30) according to claim 3 in which the tip (36) of the finger (34) is spaced radially from the center (60a) of the


axial path defined by the guide means (60) by 2.0mm to 2.5mm plus the radius of the drill bit (62) to be guided by the guide means (60).

5. A glenoid alignment guide (30) according to any of claims 1 to 4 adapted among other things for improved viewing of the glenoid face (26), the drill guide (40) having a periphery and a plurality of cut away portions (48) around the periphery to improve the view of the glenoid face (26) as the drill guide (40) is advanced. 5
6. A glenoid alignment guide (30) according to any of claims 1 to 5 in which the drill guide (40) has a generally convex-curved frontal surface (46) adapted to correspond to the generally concave-curved surface of the glenoid face (26) of the glenoid (12). 10
7. A glenoid alignment guide (30) according to any of claims 1 to 6 further comprising locking means (50) for releasably locking the drill guide (40) against the glenoid face (26) of the glenoid (12). 15
8. A glenoid alignment guide (30) according to claim 7 in which the tip (36) of the finger (34) and the guide means (60) define forward and rearward directions, respectively, along the alignment guide (30); the locking means (50) comprising: 20
  - a rod (43) projecting rearwardly from the drill guide (40);
  - a pivotable locking lever (52) having a hole (54) with an elongate cross section through which the rod (43) passes; and 25
  - resilient urging means (56) for urging the locking lever (52) such that the rod (43) is gripped by the edges of the hole (54), so that retraction of the drill guide (40) is resisted by the gripping of the rod (43) by the locking lever (52) but is possible if the locking lever (52) is manually pivoted against the urging force of the resilient urging means (56) to release the rod (43). 30
9. A glenoid alignment guide (30) according to claim 8 in which the retractor plate (32) has a track (32a) mounting the drill guide (40) on the retractor plate (32) for sliding movement of the drill guide (40) along the track (32a) in the axial direction, the track (32a) constraining the drill guide (40) from rotation relative to the retractor plate (32), the drill guide (40) being mounted on the rod (43) to move axially together with the rod (43) to allow advancing the drill guide (40) along the track (32a) of the retractor plate (32) by manually moving the rod (43). 35
10. A glenoid alignment guide (30) according to claim 9 in which the retractor plate (32) has a main body having opposite side edges, and ribs (32a) extending from the main body along the opposite side edges to define the track (32a). 40
11. A glenoid alignment guide (30) according to any of claims 1 to 10 in which the drill guide (40) has a rear surface (42), the glenoid alignment guide (30) further comprising a drill bit (62) received and guided by the drill guide (40), and a collar (66) fitted around the drill bit (62) for engaging the rear surface (42) of the drill guide (40) to limit the penetration of the drill bit (62) into the glenoid bone. 45



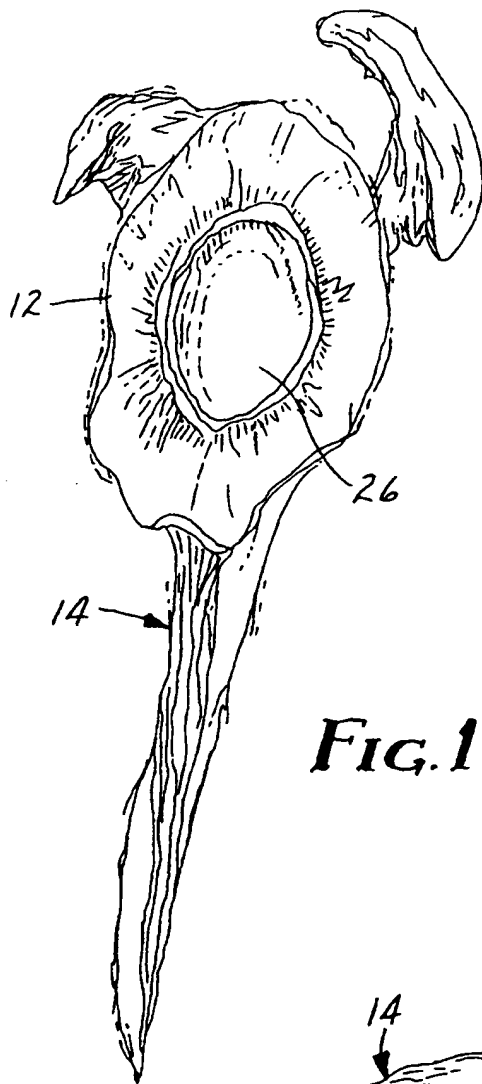

**FIG. 1**

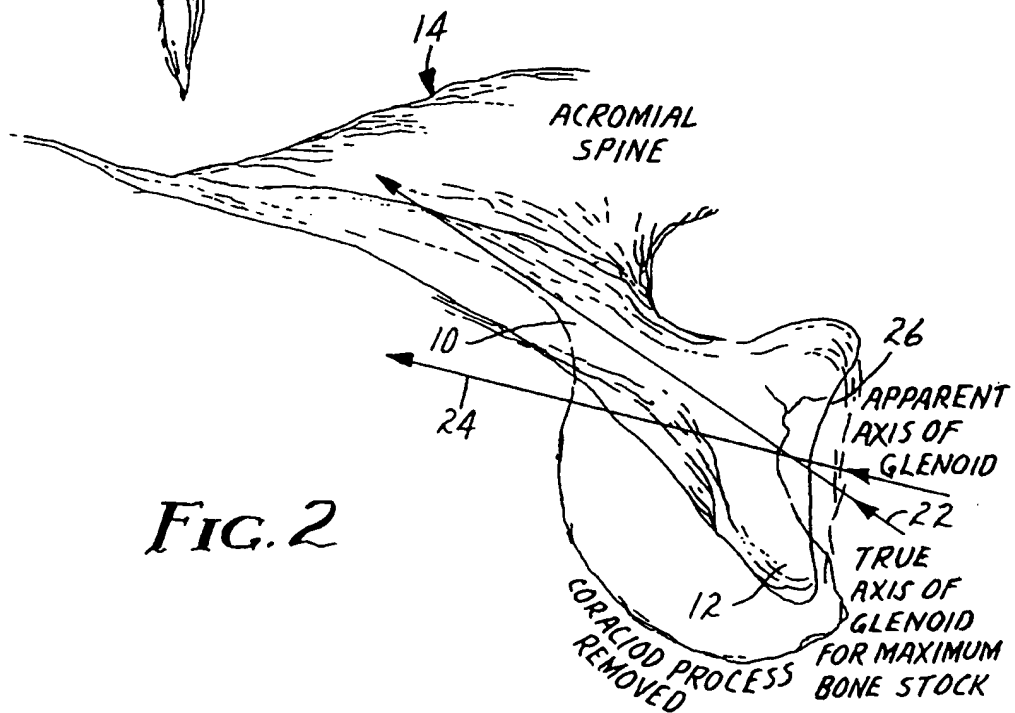

**FIG. 2**

FIG. 3





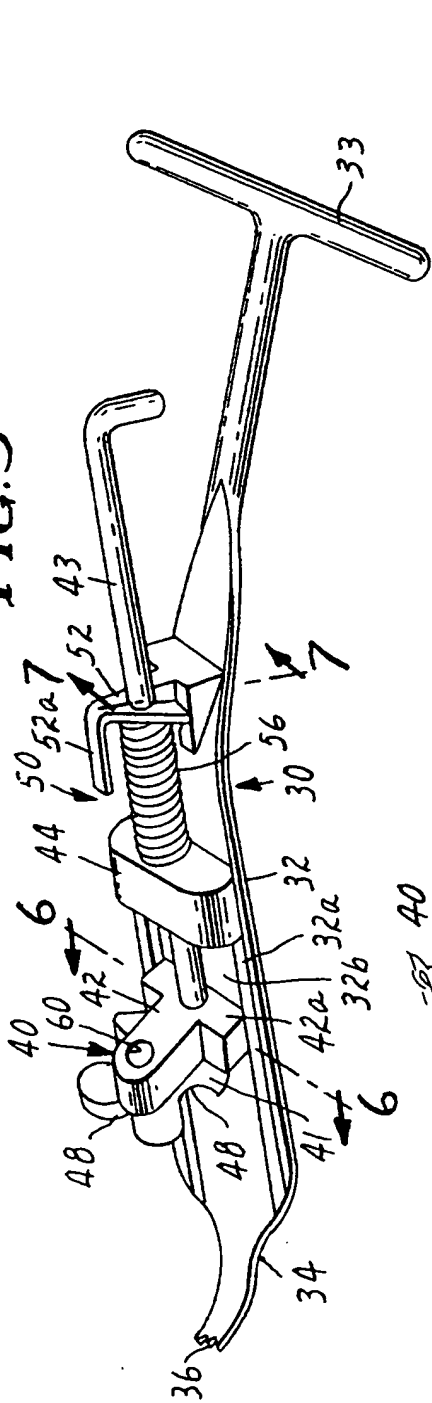


FIG. 4

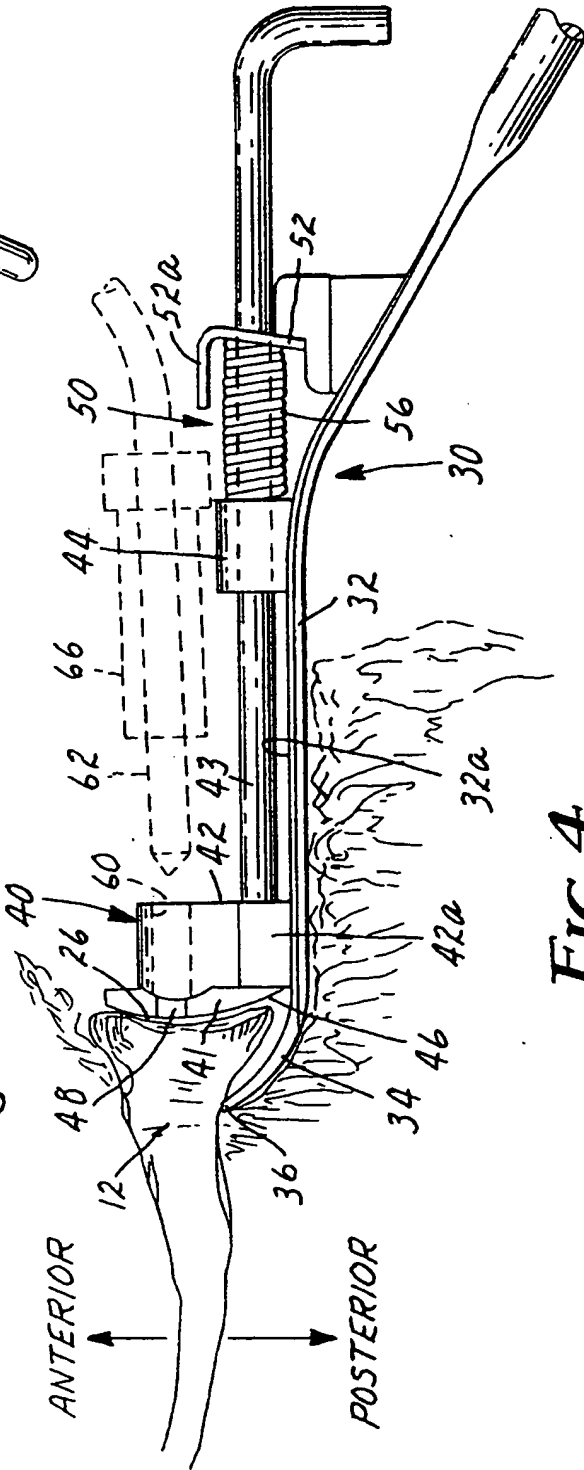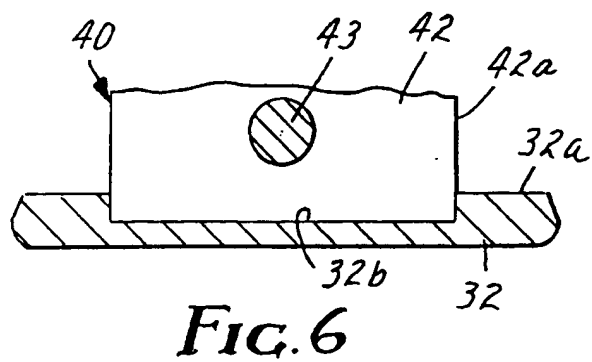

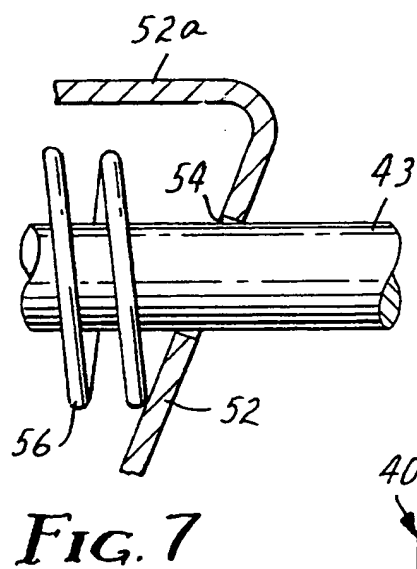

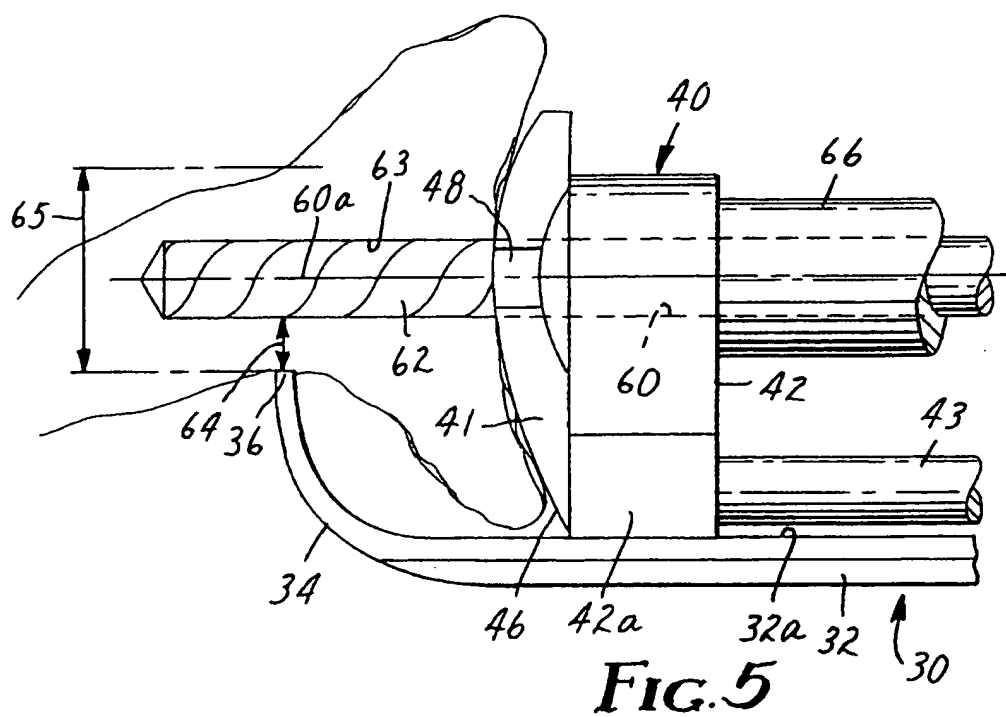

European Patent  
Office

## EUROPEAN SEARCH REPORT

Application Number  
EP 93 11 6322

DOCUMENTS CONSIDERED TO BE RELEVANT			
Category	Citation of document with indication, where appropriate, of relevant passages	Relevant to claim	CLASSIFICATION OF THE APPLICATION (Int.CLS)
A	US-A-5 030 219 (F.A.MATSEN III ET AL.) * figures 2-4 *		A61B17/16
A	--- JOURNAL OF BONE AND JOINT SURGERY. vol. 64-A, no. 2 , February 1982 , BOSTON US page 289 H.FUKUDA ET AL. 'Ring Retractor' * the whole document *		
			TECHNICAL FIELDS SEARCHED (Int.CLS)
			A61B A61F
The present search report has been drawn up for all claims			
Place of search THE HAGUE		Date of completion of the search 5 January 1994	Examiner Nice, P
<b>CATEGORY OF CITED DOCUMENTS</b> X : particularly relevant if taken alone Y : particularly relevant if combined with another document of the same category A : technological background O : non-written disclosure P : intermediate document T : theory or principle underlying the invention E : earlier patent document, but published on, or after the filing date D : document cited in the application I : document cited for other reasons & : member of the same patent family, corresponding document			